Photron



Compact high-speed cameras with high light sensitivity

The FASTCAM Mini AX is Photron's highest performance model within the FASTCAM Mini series of high-speed cameras. The Mini AX delivers exceptional light sensitivity, excellent image quality and flexible region of interest (ROI) features for customers who do not require the ultimate frame rate performance of the FASTCAM SA-X2 and SA-Z, but would benefit from the same high-end camera image sensor features.

Three performance level models - Mini AX50, AX100 and AX200 - deliver 1-megapixel image resolution (1024 x 1024 pixels) at frame rates up to 2,000fps, 4,000fps and 6,400fps respectively. All three Mini AX models offer a minimum exposure duration of 1µs as standard with recording memory options up to 32GB providing extended recording times and triggering flexibility.

Subject to export approval the Mini AX100 can be offered with maximum frame rates up to 540,000fps and the Mini AX200 with maximum frame rates up to 900,000fps with a minimum exposure time of 260 nanoseconds.

Standard operational features of the FASTCAM Mini AX include a mechanical shutter to allow remote system calibration, Gigabit Ethernet Interface for reliable system control with high-speed data transfer to PC, and the ability to remotely switch off cooling fans to eliminate vibration when recording at high magnifications.

With the combination of high frame rates, high image quality and exceptional light sensitivity contained within a 120mm x 120mm x 94mm rugged camera body weighing just 1.5kg, the FASTCAM Mini AX is ideally suited for use in a wide range of demanding scientific and industrial applications.



Model AX50 / AX100 / AX200

1-Megapixel CMOS Image Sensor: 1024 x 1024 pixels at 2,000fps (Mini AX50)

1024 x 1024 pixels at 2,000tps (Mini AX50) 1024 x 1024 pixels at 4,000fps (Mini AX100) 1024 x 1024 pixels at 6,400fps (Mini AX200)

Maximum Frame Rate:

170,000fps (Mini AX50 type 170K) 212,500fps (Mini AX100 type 200K) 540,000fps (Mini AX100 type 540K) 216,000fps (Mini AX200 type 200K) 540,000fps (Mini AX200 type 540K) 900,000fps (Mini AX200 type 900k)

Class Leading Light Sensitivity:

ISO 12232 Ssat • ISO 40,000 monochrome • ISO 16,000 color

Global Electronic Shutter:

1ms to 1µs independent of frame rate (Mini AX200 model 900K only: 260ns shutter available subject to export control)

Dynamic Range (ADC): 12-bit monochrome, 36-bit color

Compact and Lightweight:

120mm (H) x 120mm (W) x 94mm (D) 4.72" (H) x 4.72" (W) x 3.70" (D) Weight: 1.5Kg (3.30 lbs.)

Internal Recording Memory: 8GB, 16GB, or 32GB

Fast Gigabit Ethernet Interface:

Provides high-speed image download to standard notebook/PC

Flexible Frame Synchronization:


Frame rate may be synchronized to external unstable frequencies

High-G Rated:

Suitable for application in high-G environments; operation tested to 100G, 10ms, 6-axes

Fan Stop Function:

Remotely switch off cooling fans to eliminate vibration

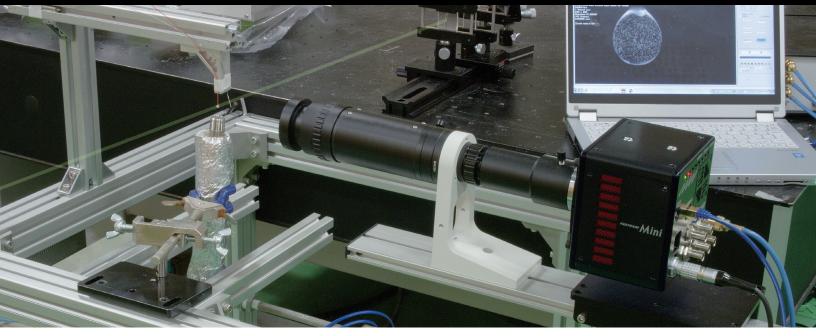

Light Sensitivity:

Expressions of light sensitivity in high-speed cameras can be confusing as a variety of differing measurement techniques are used. Photron publishes light sensitivity figures for its products using the ISO 12232 Ssat Standard.

| FASTCAM MINI AX | ISO 12232 Ssat |
|-------------------|----------------|
| Monochrome models | ISO 40,000 |
| Color models | ISO 16,000 |

ISO 12232 Ssat values published by Photron for both monochrome and color cameras are measured excluding infrared sensitivity as defined by the ISO standard measurement procedure ISO 14524.

Monochrome sensors used in the FASTCAM Mini AX cameras are supplied without an IR absorbing filter, extending the camera spectral response beyod 900nm. When the sensitivity of the FASTCAM Mini AX camera is measured to tungsten light including near IR response an equivalent value of ISO 100,000 is obtained.

Image Sensor:

The FASTCAM Mini AX system uses an advanced CMOS image sensor optimized for light sensitivity and high image quality that is unique to Photron.


A 20-micron pixel pitch gives a sensor size at full image resolution of 20.48 x 20.48mm (diagonal 28.96mm).

Lenses designed for both FX (35mm full frame) and also DX (APS-C digital SLR) formats are fully compatible with the FASTCAM Mini AX at full image resolution.

| Sensor Type | Proprietary Design Advanced CMOS | | |
|-----------------------------|---|--|--|
| Maximum Resolution (pixels) | 1024 x 1024 pixels | | |
| Sensor Size / Diagonal | 20.48 x 20.48mm / 28.96mm | | |
| Pixel Size (microns) | 20µm x 20µm | | |
| Quantum Efficiency | 46% at 630nm | | |
| Fill Factor | 58% | | |
| Color Matrix | Bayer CFA (single sensor) | | |
| ISO 12232 Ssat sensitivity | ISO 40,000 monochrome ISO 16,000 color (monochrome sensor equivalent ISO 100,000 including near IR response) | | |
| Shutter | Global Electronic Shutter 1ms to 1µs independent of frame rate (Mini AX200 model 900K only; 260ps shutter | | |

Shutter

frame rate (Mini AX200 model 900K only: 260ns shutter available subject to export control)

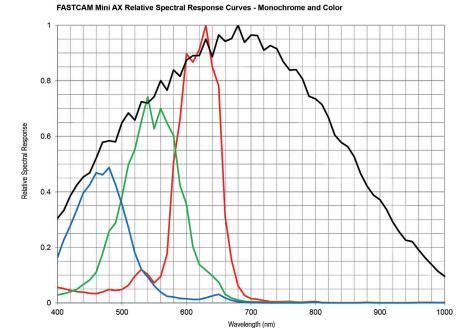

Camera Performance Specifications

| Model | Mini AX50 | Mini AX100 | Mini AX200 | |
|--|--|---|---|--|
| Full Frame Performance | 2,000fps 1024 x 1024 pixels | 4,000fps 1024 x 1024 pixels | 6,400fps 1024 x 1024 pixels | |
| Maximum Frame Rate | Type 170K: 170,000fps (128 x 16 pixels) | Type 200K: 212,500fps (128 x 16 pixels) Type 540K: 540,000fps* (128 x 16 pixels) | Type 200K: 216,000fps (128 x 16 pixels) Type 540K: 540,000fps* (128 x 16 pixels) Type 900K: 900,000fps* (128 x 16 pixels) | |
| Minimum Exposure Time | Global electronic shutter to 1.05µs selectable | e independent or frame rate (260ns option availa | ble with Mini AX200 type 900K only) * | |
| Inter Frame Time (for PIV) | 1.71µs | | | |
| Ruggedized Mechanical Calibration Shutter | Standard feature | | | |
| Dynamic Range (ADC) | 12-bit monochrome 36-bit color | | | |
| Memory Capacity Options | 8GB: 5,457 frames at full resolution 16GB: 10,918 frames at full resolution 32GB: 21,841 frames at full resolution | | | |
| Memory Partitions | Up to 64 memory segments | | | |
| Region of Interest | Selectable in steps of 128 pixels (horizontal) x 16 pixels (vertical) | | | |
| Trigger Inputs | Selectable +/- TTL 5V and switch closure | | | |
| Trigger Delay | Programmable on selected input / output trig | gers: 100ns resolution | | |
| Input / Output | Input: Trigger (TTL/Switch), sync, ready, eve Output: trigger, sync, ready, rec, exposure | nt, IRIG | | |
| Trigger Modes | Start, end, center, manual, random, random reset, image trigger, time lapse | | | |
| Time Code Input | IRIG-B | | | |
| External Sync | +/- TTL 5Vp-p Variable frequency sync | | | |
| Camera Control Interface | High-speed Gigabit Ethernet | | | |
| Image Data Display | Frame rate, shutter speed, trigger mode, date/time, status, real time / IRIG time, frame count, resolution | | | |
| Saved Image Formats | BMP, TIFF, JPEG, PNG, RAW, RAWW, MRAW, AVI, WMV, FTIF, MOV - Images can be saved with or without image data and in 8-bit, 16-bit or 36-bit depth of sensor where supported | | | |
| Supported OS | Microsoft Windows operating system including: 7, 8, 8.1, 10 (32/64-bit) | | | |

* Frame rates above 225,000 fps and exposure times below 1µs may be subject to export control regulations in some areas

High-Speed Gigabit Ethernet Interface:

The FASTCAM Mini AX camera system is equipped with a high-speed Gigabit Ethernet Interface to provide reliable network communication and fast download of image data.

Dedicated I/O:

A dedicated BNC connection for a contact closure hardware trigger input is provided. In addition, two programmable inputs and two programmable output channels provide direct connection for common tasks such as synchronization of multiple cameras and operation in conjunction with Data Acquisition (DAQ) hardware.

High-G Mechanical Calibration Shutter:

The ruggedized mechanical shutter fitted as standard to the FASTCAM Mini AX camera allows sensor black balance calibration to be carried out remotely from the system control software.

Nikon G-Type Compatible Lens Fitting:

The FASTCAM Mini AX camera is equipped with an objective lens mount compatible with readily available Nikon G-type lenses. Controls provided within the lens mount allow the control of lens aperture on lenses without external iris control.





Camera Operation Features

| Frame Synchronization | Accurate frame synchronization with other cameras and with external and unstable frequencies. | |
|--|--|--|
| Dual Slope Shutter (Extended Dynamic Range) | Selectable in 20 steps (0 to 95% in 5% increments) to prevent pixel overexposure without post processing. | |
| Memory Partitions | Up to 64 memory segments allow multiple events to be stored in camera memory before downloading, with automatic progression to the next available partition. | |
| Low Light Mode | Operation at minimum frame rate with separately adjustable shutter time to allow easy camera set-up and focus in ambient lighting. | |
| IRIG Phase Lock | Enables multiple cameras to be synchronized together with other instrumentation equipment or to a master external time source. | |
| Internal Time Delay Generator | Allows programmable delays to be set on input and output triggers; 100ns resolution. | |
| Event Markers | Up to ten user-entered event markers to define specific events within the recorded image sequence . | |
| Download While Recording | FASTCAM Mini AX supports Partition Recording Mode, allowing image data captured in one memory partition to be downloaded while at the same time recording into another partition. | |
| Automatic Download | The system can be set to automatically download image data to the control PC and, when download is complete to re-arm in readiness for the next trigger with automatically incremented file names. | |
| Software Binning | Virtual pixel binning (2x2, 4x4 etc.) allows increased light sensitivity with reduced image resolution without changing camera field of view. | |

Operation Software Features

| Image Calibration | 2D image calibration allows the measurement of distance and angle from the image. A calibration grid overlay can be superimposed on the image. |
|---|---|
| Image Overlay | A stored reference image may be overlaid on the live image to allow accurate camera positioning to achieve the same view as a previous test. |
| Import of Multiple Image Sequences | Multiple image sequences can be loaded and simultaneously replayed. Timing of image sequences can be adjusted to create a common time reference. Time based synchronization allows images captured at different frame rates to be synchronized. |
| High Dynamic Range Mode Making use of the full sensor dynamic range, HDR mode allows enhanced detail in both light and dark areas of an image to be a simultaneously. | |
| Motion Detector In order to highlight subtle changes in an image, Motion Detector allows a reference image to be subtracted from a recorder including propagation of shock waves and surface changes during impact can be visualized using the feature. | |
| Line Profile A line profile representing grey levels along a line drawn across any region of the image is displayed. In live mode the Line Profile ensure optimum image focus is achieved. | |
| Histogram | A histogram displaying grey levels within a user-defined image area is displayed. In live mode the Histogram can be used to ensure that optimum exposure levels are set for the scene being recorded. |

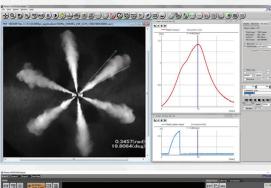

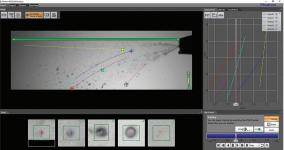

Photron FASTCAM Viewer:

Photron FASTCAM Viewer software (PFV) has been designed to provide an intuitive and feature rich user interface for the control of Photron high-speed cameras, data saving, image replay and simple motion analysis. Advanced operation menus provide access to features for advanced camera operation and image enhancement. Tools are provided to allow image calibration and easy measurement of angles and distances from image data. Also included are a C++ SDK and wrappers for LabView and MATLAB ®.

An optional software plug-in module provides synchronisation between Photron high- speed cameras and data acquired through National Instruments data acquisition systems. Synchronised data captured by the DAQ system provides waveform information which can be viewed alongside high-speed camera images.

Photron FASTCAM Analysis:

PFV software allows image sequences to be exported directly to optional Photron FASTCAM Analysis (PFA) Motion Analysis software. This entry level Motion Analysis software with an on screen 'step by step guide' function launches automatically from Photron FASTCAM Viewer software, and provides automated tracking of up to 5 points using feature or correlation tracking algorithms for the automated analysis of motion within an image sequence.

Variable Region of Interest:

Region of Interest (ROI) or sub-windowing allows a user-specified portion of the sensor to be defined to capture images. By using a reduced portion of the image area, the frame rate at which images are recorded can be increased. FASTCAM Mini AX allows the ROI to be set in increments of 128 pixels horizontal and 16 pixels vertical.

Square Image Sensor Format:

Unlike broadcast and media applications where image formats such as 16:9 have now become standard, in scientific and industrial imaging applications an image sensor with a 1:1 image format is generally accepted to be advantageous. To capture the maximum useful image data in applications including microscopy, detonics, combustion imaging and many others, a 1:1 sensor format provides greater flexibility than 'letterbox' image formats. The FASTCAM Mini AX image sensor allows the user to choose either square or rectangular image formats in order to obtain the maximum subject information.

External Frame Synchronization:

The FASTCAM Mini AX camera can be fully synchronized with an external event to allow the timing of when each individual image is captured to be precisely referenced. The camera can be accurately synchronized to unstable frequencies allowing complex events such as combustion in rapidly accelerating or decelerating engines to be recorded and studied.

Record During Download Operation:

FASTCAM Mini AX recording memory can be divided into multiple active sections. The user can record an on-going event in one memory partition while at the same time downloading a previously recorded image sequence in order to improve workflow and optimize camera operation.



| | Mini AX200 | | | | | | |
|----------------|------------|-----------|--------------|-----------|--------------|------------|--------------|
| Resolution | Frame Rate | 8GB | | 16GB | | 32GB | |
| (h x v pixels) | Max fps | Frames | Time (sec)** | Frames | Time (sec)** | Frames | Time (sec)** |
| 1024 x 1024 | 6,400 | 5,457 | 0.85 | 10,918 | 1.71 | 21,841 | 3.41 |
| 1024 x 896 | 7,200 | 6,236 | 0.87 | 12,478 | 1.73 | 24,961 | 3.47 |
| 896 x 896 | 8,100 | 7,127 | 0.88 | 14,261 | 1.76 | 28,527 | 3.52 |
| 768 x 768 | 10,800 | 9,701 | 0.90 | 19,410 | 1.80 | 38,829 | 3.60 |
| 512 x 512 | 22,500 | 21,829 | 0.97 | 43,674 | 1.94 | 87,365 | 3.88 |
| 512 x 256 | 43,200 | 43,658 | 1.01 | 87,349 | 2.02 | 174,730 | 4.04 |
| 256 x 256 | 67,500 | 87,317 | 1.29 | 174,698 | 2.59 | 349,461 | 5.18 |
| 256 x 128 | 120,000 | 174,634 | 1.46 | 349,397 | 2.91 | 698,922 | 5.82 |
| 128 x 128 | 162,000 | 349,269 | 2.16 | 698,794 | 4.31 | 1,397,845 | 8.63 |
| 128 x 64 | 259,200 | 698,538 | 2.69 | 1,397,589 | 5.39 | 2,795,690 | 10.79 |
| 128 x 32 | 360,000 | 1,397,077 | 3.88 | 2,795,178 | 7.76 | 5,591,381 | 15.53 |
| 128 x 16 | 540,000 | 2,794,154 | 5.17 | 5,590,357 | 10.35 | 11,182,762 | 20.71 |
| 128 x 16 | 900,000 | | | | | | |

Mini AX100


| Resolution | Frame Rate | 8GB | | 16GB | | 32GB | |
|----------------|------------|-----------|--------------|-----------|--------------|------------|--------------|
| (h x v pixels) | Max fps | Frames | Time (sec)** | Frames | Time (sec)** | Frames | Time (sec)** |
| 1024 x 1024 | 4,000 | 5,457 | 1.36 | 10,918 | 2.73 | 21,841 | 5.46 |
| 1024 x 896 | 4,500 | 6,236 | 1.39 | 12,478 | 2.77 | 24,961 | 5.55 |
| 896 x 896 | 5,400 | 7,127 | 1.32 | 14,261 | 2.64 | 28,527 | 5.28 |
| 768 x 768 | 6,800 | 9,701 | 1.43 | 19,410 | 2.85 | 38,829 | 5.71 |
| 512 x 512 | 13,600 | 21,829 | 1.61 | 43,674 | 3.21 | 87,365 | 6.42 |
| 512 x 256 | 25,500 | 43,658 | 1.71 | 87,349 | 3.43 | 174,730 | 6.85 |
| 256 x 256 | 37,500 | 87,317 | 2.33 | 174,698 | 4.66 | 349,461 | 9.32 |
| 256 x 128 | 61,200 | 174,634 | 2.85 | 349,397 | 5.71 | 698,922 | 11.42 |
| 128 x 128 | 76,500 | 349,269 | 4.57 | 698,794 | 9.13 | 1,397,845 | 18.27 |
| 128 x 64 | 127,500 | 698,538 | 5.48 | 1,397,589 | 10.96 | 2,795,690 | 21.93 |
| 128 x 32 | 170,000 | 1,397,077 | 8.22 | 2,795,178 | 16.44 | 5,591,381 | 32.89 |
| 128 x 16 | 540,000 | 2,794,154 | 5.17 | 5,590,357 | 10.35 | 11,182,762 | 20.71 |

Mini AX50 16GB Resolution Frame Rate 8GB 32GB Frames (h x v pixels) Max fps Frames Time (sec)** Time (sec)* Frames Time (sec)** 1024 x 1024 2,000 5,457 10,918 21,841 2.73 5.46 10.92 1024 x 896 2,500 6,236 2.49 12,478 4.99 24,961 9.98 2,500 896 x 896 7,127 2.85 14,261 5.70 28,527 11.41 768 x 768 3,600 9,701 2.69 19,410 38,829 10.79 5.39 512 x 512 7,200 21,829 43,674 3.03 6.07 87,365 12.13 512 x 256 13,600 43,658 3.21 87,349 6.42 174,730 12.85 256 x 256 20,400 87,317 4.28 174,698 8.56 349,461 17.13 256 x 128 37,500 174,634 4.66 349,397 9.32 698,922 18.64 128 x 128 45,900 349,269 7.61 698,794 15.22 30.45 1,397,845 9.13 36.54 128 x 64 76,500 698,538 1,397,589 18.27 2,795,690 128 x 32 127,500 1,397,077 10.96 2,795,178 21.92 5,591,381 43.85 128 x 16 170,000 2,794,154 16.44 5,590,357 32.88 11,182,762 65.78

* Specifications subject to change without notice.

** Recording time is an estimate and may be different depending on recording conditions and settings.

Photo Schlieren imaging of fuel injection and engine combustion 20,000fps

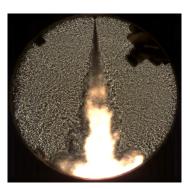

Compatibility with Specialist Lens Systems:

A combination of small physical size, low weight and high light sensitivity allows the FASTCAM Mini AX to be coupled to a range of optical systems such as scientific and long distance microscopes, rigid endoscopes or borescopes, and image intensifiers for applications ranging from imaging flows in microfluidic devices to combustion diagnostics.

PIV and DIC Requirements:

FASTCAM Mini AX specifications match with the requirements for optical measurement techniques such as Particle Image Velocimetry (PIV) and Digital Image Correlation (DIC). The FASTCAM Mini AX has many key performance specifications desired for these measurement systems.


In PIV the detection of low light levels from small particles is fundamental. A high sensitivity image sensor allows the use of smaller tracer particles and/ or lower laser power.

For DIC applications a highly sensitive camera allows the use of smaller objective lens apertures yielding greater depth of field and enhanced measurement of out of plane displacements.

Small Physical Size:

The small physical size and weight of the Mini camera range allows the use of conventional opto-mechanical hardware for rigid and stable mounting of multiple cameras, and for the location of cameras in space limited locations.

Specifications subject to change without notice.

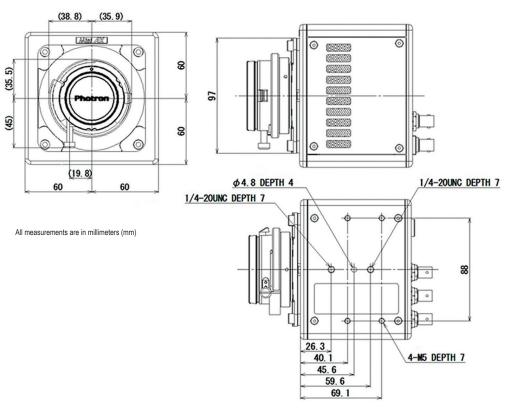

Mechanical and Environmental Specifications

| Mechanical | |
|--|---|
| Lens Mount | F-mount (G-type lens compatible) and C-mount provided - Optional lens mounts available include M42 adapter |
| Camera Mountings | 4 x 1/4 - 20 UNC (base and top), 4 x M5 (base) |
| External Dimensions | |
| Camera Body (excluding protrusions) | 120mm (H) x 120mm (W) x 94mm (D) 4.72" (H) x 4.72" (W) x 3.70" (D) |
| Weight | |
| Camera Body | 1.5kg (3.30lbs) |
| Environmental | |
| Operating Temperature | 0 to 40C, 32° to 104°F |
| Storage Temperature | -20 to 60C, -4° to 140°F |
| Humidity | 85% or less (non-condensing) |
| Cooling | Internal fan cooling (fan-off mode supported) |
| Operational Shock | 100G, 10ms, 6-axes |
| Power | |
| AC Power (with supplied adapter) | 100 to 240V, 50 to 60Hz |
| DC Power | 22 to 32V, 55VA |

PHOTRON USA, INC. 9520 Padgett Street, Suite 110 San Diego, CA 92126 USA

Tel: 858.684.3555 or 800.585.2129 Fax: 858.684.3558 Email: image@photron.com www.photron.com PHOTRON EUROPE LIMITED The Barn, Bottom Road West Wycombe, Bucks. HP14 4BS United Kingdom

Tel: +44 (0) 1494 481011 Fax: +44 (0) 1494 487011 Email: image@photron.com www.photron.com PHOTRON (Shanghai) Room 20C, Zhao-Feng World Trade Building No. 369, JiangSu Road ChangNing District Shanghai, 200050 China Tel: +86 (0) 21-5268-3700 Email: info@photron.cn.com

PHOTRON LIMITED

Kanda Jinbo-cho 1-105 Chiyoda-ku, Tokyo 101-0051 Japan

Tel: +81 (0) 3 3518-6271 Fax: +81 (0) 3 3518-6279 Email: image@photron.co.jp www.photron.co.jp